



Oval, POLYsmooth™



Oval, MESMO®



Round, POLYsmooth™

Gluteus Implants

Silicone-gel-filled gluteus implants available in two shapes and surfaces, with planar or concave back

Implants made by POLYTECH – QUALITY made in Germany



Gluteus Implants, € 0483

Product Description

Silicone gel-filled gluteus implants for compensation of soft tissue deficiencies are available with two different surfaces:

- POLYsmooth™ (smooth)
- MESMO® (micro-textured)

The surgeon can choose among oval and round shapes of various sizes as well as between a concave or planar implant posterior. Thus he is in a position to achieve the optimal result for either augmentation or reconstruction of the gluteus region in patients of both sexes.

Information on Materials

The shell of all gluteus implants consists of a strong, resilient and highly elastic silicone elastomer equipped with a barrier layer. It is filled with a highly cross-linked silicone gel. Due to the high stress that can possibly act on the implants, a gel with firm consistency has been chosen. The surface of the implants is either smooth, POLYsmooth™, or micro-textured, MESMO®. The MESMO® texture features a small pore size of an average 100 to 150µm.

Production of the implants is based on a quality management according to ISO 9001 and ISO 13485. All materials used are of implantable medical quality.

Indications

- Augmentation
- Contour correction
- Reconstruction, e.g. after trauma
- Implant exchange

Important Information

Gluteal implants are delivered individually packed and sterile.


Attention: Please note the detailed instructions for use included with each product.

Gluteus implants – concave posterior

Silicone-gel filled, individually packed, sterile, € 0483

| Order n° | Vol. (ml) | Dimensions (mm) | | | Posterior |
|----------|-----------|-----------------|---|---|-----------|
| | | A | B | C | |

| Oval, POLYsmooth™ | | | | | |
|------------------------|-----|-----|-----|----|---------|
| 10300-115 ¹ | 115 | 80 | 120 | 27 | concave |
| 10300-215 ¹ | 215 | 105 | 142 | 30 | concave |
| 10300-290 ¹ | 290 | 114 | 152 | 33 | concave |
| 10300-360 ¹ | 360 | 125 | 163 | 34 | concave |
| 10300-450 ¹ | 450 | 132 | 174 | 38 | concave |

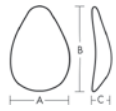

| Oval, MESMO® | | | | | |
|------------------------|-----|-----|-----|----|---------|
| 15300-115 ¹ | 115 | 80 | 120 | 27 | concave |
| 15300-215 ¹ | 215 | 105 | 142 | 30 | concave |
| 15300-290 ¹ | 290 | 114 | 152 | 33 | concave |
| 15300-360 ¹ | 360 | 125 | 163 | 34 | concave |
| 15300-450 ¹ | 450 | 132 | 174 | 38 | concave |



| Round, POLYsmooth™ | | | | | |
|------------------------|-----|-----|-----|----|---------|
| 10301-200 ¹ | 200 | 107 | 107 | 35 | concave |
| 10301-250 ¹ | 250 | 112 | 112 | 41 | concave |
| 10301-300 ¹ | 300 | 118 | 118 | 44 | concave |
| 10301-360 ¹ | 360 | 127 | 127 | 43 | concave |




Gluteus implants – planar posterior


Silicone-gel filled, individually packed, sterile, € 0483

| Order n° | Vol. (ml) | Dimensions (mm) | | | Posterior |
|----------|-----------|-----------------|---|---|-----------|
| | | A | B | C | |

| Oval, POLYsmooth™ | | | | | |
|------------------------|-----|-----|-----|----|--------|
| 10300-240 ¹ | 240 | 90 | 130 | 40 | planar |
| 10300-305 ¹ | 305 | 95 | 139 | 45 | planar |
| 10300-365 ¹ | 365 | 100 | 148 | 48 | planar |
| 10300-420 ¹ | 420 | 109 | 152 | 50 | planar |

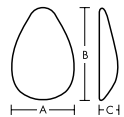

| Oval, MESMO® | | | | | |
|------------------------|-----|-----|-----|----|--------|
| 15300-240 ¹ | 240 | 90 | 130 | 40 | planar |
| 15300-305 ¹ | 305 | 95 | 139 | 45 | planar |
| 15300-365 ¹ | 365 | 100 | 148 | 48 | planar |
| 15300-420 ¹ | 420 | 109 | 152 | 50 | planar |

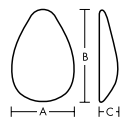

| Round, POLYsmooth™ | | | | | |
|------------------------|-----|-----|-----|----|--------|
| 10301-225 ¹ | 225 | 107 | 107 | 35 | planar |
| 10301-275 | 275 | 112 | 112 | 41 | planar |
| 10301-340 | 340 | 118 | 118 | 44 | planar |
| 10301-390 ¹ | 390 | 128 | 128 | 44 | planar |
| 10301-560 ¹ | 560 | 141 | 141 | 49 | planar |



approximate measurements

| | |
|----------|---|
| On stock | ¹ Supply after order ≈ 6 weeks |
|----------|---|